

HIRU Start Guide

Mounting system - Holder

Hiru | The first multi-platform
eye tracker in the world.




Holder

Hiru can be used with a [laptop](#) or [screen](#). For these cases, Hiru can be used with our holder.



If you will be using a holder, you will receive the following [elements](#) in your box: 2 [screws](#), 2 [magnetic strips](#) and the [holder](#).

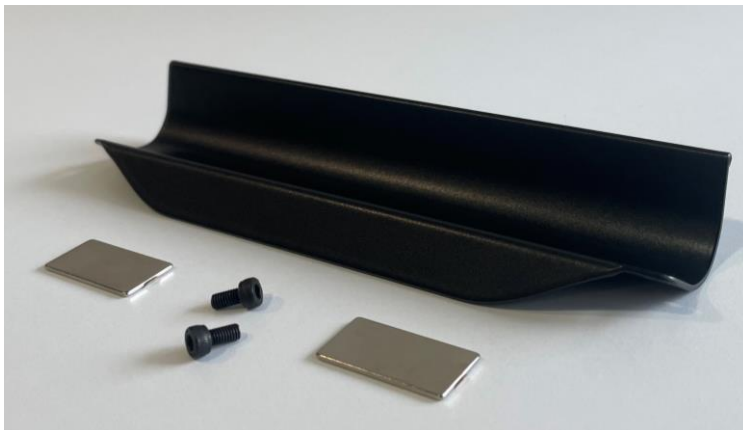

Holder

Using the Hiru with a laptop

Screw the Hiru to the holder. The screwed threads of the Hiru are located in the bottom side.



Once screwed, the Hiru will be correctly located in the holder and the working angles will be optimum. Locate the Hiru with the holder on the laptop, below the screen.



Holder

Using the Hiru with a screen

For using Hiru with a screen, the Hiru and the holder will be attached to the screen using the magnetic strips.

For locating the magnetic strips correctly, follow the next steps:

1. Locate the magnetic strips on the holder, below the supporting fold.
2. Clean the area of the monitor where the magnetic strips will be adhered.
3. Take away the protection strip of the adhesives and locate them in the lower side of the monitor. Keep pressing for 2 minutes.
4. Let the adhesive attach correctly for 1-2 hours.
5. After the resting time, remove the holder and screw the Hiru in the holder for getting the optimum angle.
6. If you want to use Hiru in another situation, just remove the holder and Hiru from the screen, the magnetic strips will remain there for next occasions.

